Paths of analysis*

Synthia

March 3, 2022

Analysis parameters 1

Analysis type: Automatic Retrosynthesis

Rules: none selected

Filters: FGI, FGI with protections

Max. paths returned: 5

Max. iterations: 300

Commercial:

1. Max. molecular weight - 1000 g/mol

2. Max. price - 1000 \$/g

Published:

- 1. Max. molecular weight 1000 g/mol
- 2. Popularity 10

My Stockroom:

1. Max. molecular weight - 1000 g/mol

Reaction scoring formula: TUNNEL COEF*FGI COEF*STEP*20+1000 000*(CONFLICT+NON SELECTIVITY+FILTERS+PROTECT)

Chemical scoring formula: SMALLER^ 3,SMALLER^ 1.5

Min. search width: 400

Max. reactions per product: 60

Strategies: none selected

^{*}The results stated herein were generated using the proprietary platform owned and maintained by Grzybowski Scientific Inventions, Inc., a subsidiary of Merck KGaA, Darmstadt Germany. The results are provided on an as is basis, and shall be used solely in connection with the rights afforded in the license agreement and for no other purpose.

FGI Coeff: 0

JSON Parameters: {}

2 Paths

2 paths found. Paths are sorted by score. Reactions are sorted in appearance order for each path.

2.1 Path 1

Score: 101.85


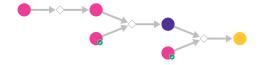


Figure 1: Outline of path 1

2.1.1 Balz-Schiemann Reaction

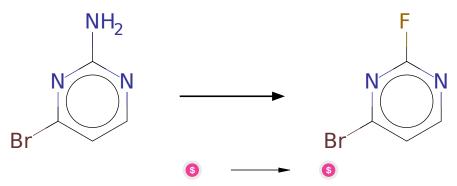

Substrates:

1. 4-Bromopyrimidin-2-amine - Combi-Blocks

Products:

1. 4-Bromo-2-fluoropyrimidine - Enamine

Typical conditions: NaNO2.HF-pyridine.-25 to 0C

Protections: none
Yield: moderate

Reference: 10.1021/jm100432w and 10.1021/jm5008177 and 10.1021/o1401540k and 10.1021/jm401551n (main text and SI, page S8) and 10.1515/chempap-2016-0033 and 10.1021/jm0311442 and 10.1016/j.jfluchem.2007.03.012 and 10.1021/jo00185a023

Retrosynthesis ID: 29906

2.1.2 Meerwein coupling of diazonium salt with heteroaryl

Substrates:

1. 4-Bromo-2-fluoropyrimidine - Enamine

2. 3-Chloroaniline - available at Sigma-Aldrich

Products:

1. Fc1ncc(Br)c(-c2ccc(Cl)c2)n1

Typical conditions: 1) HCl.NaNO2 2) [Ru(bpy)3Cl2]*6H2O.45W

bulb. H2O.rt

Protections: none
Yield: moderate

Reference: 10.1002/chem.201304120

2.1.3 Nucleophilic aromatic substitution

$$H_{2}N$$
 $+$
 CI
 Br
 CI
 Br
 CI
 Br

Substrates:

- 1. Fc1ncc(Br)c(-c2cccc(Cl)c2)n1
- 2. 5-((4-Ethylpiperazin-1-yl)methyl)pyridin-2-amine $available\ at\ Sigma-Aldrich$

Products:

 $1. \ \ CCN1CCN(Cc2ccc(Nc3ncc(Br)c(-c4cccc(Cl)c4)n3)nc2)CC1$

Typical conditions: Solvent

Protections: none

Yield: good

Reference: 10.1002/9781118093559.ch4

Retrosynthesis ID: 49476

2.2 Path 2

Score: 119.28

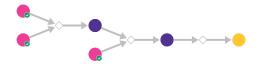


Figure 2: Outline of path 2

${\bf 2.2.1} \quad {\bf Substitution \ of \ 2-thiomethyl pyrimidines \ with \ amines}$

Substrates:

- 1. 4-Bromo-2-(methylthio)pyrimidine available at Sigma-Aldrich
- 2. 5-((4-Ethylpiperazin-1-yl)methyl)pyridin-2-amine available at Sigma-Aldrich

Products:

 $1. \ CCN1CCN(Cc2ccc(Nc3nccc(Br)n3)nc2)CC1 \\$

Typical conditions: K2CO3.DMF

Protections: none

Yield: good

Reference: 10.1021/jm980222w AND 10.1016/j.cclet.2014.10.007 AND

 $10.1002/jhet.5570280520 \ AND \ 10.1080/00397910701396930$

2.2.2 Suzuki coupling of arylboronic acids with aryl bromides

Substrates:

1. CCN1CCN(Cc2ccc(Nc3nccc(Br)n3)nc2)CC1

2. 3-Chlorophenylboronic acid - available at Sigma-Aldrich

Products:

1. CCN1CCN(Cc2ccc(Nc3nccc(-c4cccc(Cl)c4)n3)nc2)CC1

 ${\bf Typical\ conditions:}\ {\bf Pd\ catalyst.base.solvent}$

Protections: none

Yield: good

Reference: 10.1021/cr00039a007 and $10.1007/3418_2012_32$ and 10.1021/cr0505268 and 10.1016/j.jfluchem.2016.01.018 and 10.1039/C3CS60197H and 10.1016/j.ejmech.2018.08.092 and 10.1038/s41929-020-00564-z (metal-free coupling)

2.2.3 Bromination of aromatic compounds

Substrates:

 $1. \ \ CCN1CCN(Cc2ccc(Nc3nccc(-c4cccc(Cl)c4)n3)nc2)CC1$

Products:

 $1. \ CCN1CCN(Cc2ccc(Nc3ncc(Br)c(-c4cccc(Cl)c4)n3)nc2)CC1 \\$

Typical conditions: Br2.Fe

Protections: none

Yield: good

Reference: 10.1021/acs.accounts.6b00120